



DW-HL528S

M-Cascade system ULT freezer that consists of two stages cascade cooling technology with absolute hydrocarbon refrigerant, which are R290 and R170 as natural gas. The inverter compressors programmed to help energy saving under long period preservation or provide fast temperature recovery while frequently door open to maintain specimens liveness.



Energy consumption



Pulling down to -80°C



Temp. recovery after door opened for 1 min

*Typical freezer data based on internal testing with freezer set point at -80°C and ambient temperature at 25°C.

Features



Cascade system
Two stages cascade system



Touch screen
10"LCD touch screen



Hydrocarbon refrigerant
Environment-friendly natural gas



Inverter technology
Two inverter compressors




Vacuum insulation
VIP insulation foaming



Alarms system
12 kinds of alarm functions



Structure

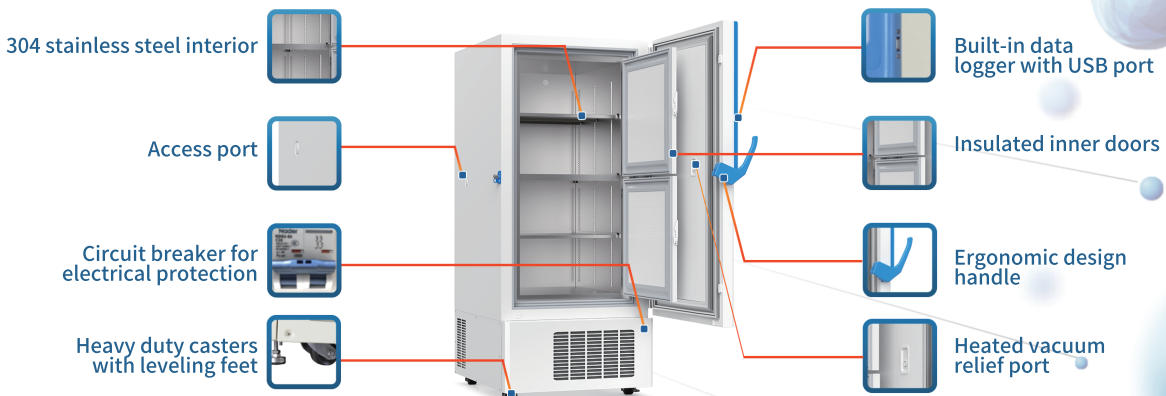



Dimensions

Specification

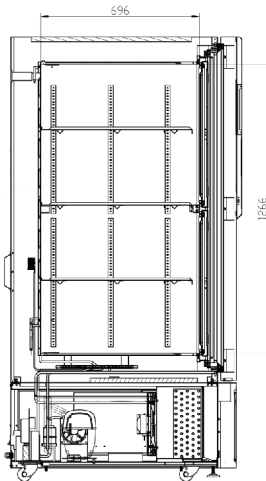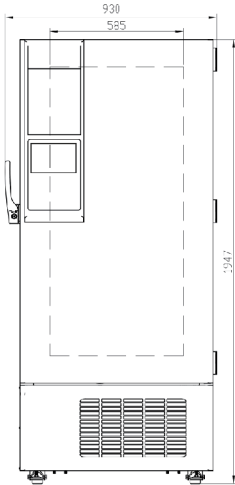

-86°C Ultra-low Temperature Freezer	
Model	DW-HL528S
Capacity	528L/18.7Cu.Ft.
Internal Size (W*D*H) mm/inch	585*696*1266/23.0*27.4*49.8
External Size (W*D*H) mm/inch	930*1041*1947/36.6*41.0*76.7
Package Size (W*D*H) mm/inch	1020*1166*2132/40.16*45.91*83.94
NW/GW (kgs)	260/310
Performance	
Temperature Range	-40~-86°C
Ambient Temperature	16-32°C
Cooling Performance	-86°C
Climate Class	N
Controller	Microprocessor
Display	10" HD Intelligent Touch Screen
Refrigeration	
Compressors	2pcs
Cooling Method	Direct Cooling
Defrost Mode	Manual
Refrigerants	High Stage: R290, Low Stage: R170
Insulation Thickness(mm)	D:117,L/R:130,B:130
Construction	
External Material	Cold Rolled Steel Sheet
Inner Material	Stainless Steel
Shelves	3(stainless steel)
Access Test Port	2pcs. Ø 25 mm
Casters	4+(2 balancing feet)
Data Logging/ Interval /Recording Time	USB/Record every 5 minutes / 10 years
Interface	Remote Alarm Contact, RS485
Backup Battery	Yes
Alarms	
Temperature	High/Low Temperature, High Ambient Temperature
Electrical	Power Failure, Low Battery, High-Stage High Temperature Alarm, Battery Level Detection Fault Alarm
System	Sensor Error, Main Board Communication Error, Door Ajar, High Pressure Alarm, Condenser Overheating Alarm
Accessories	
Freezer Rack Model	ZKL304-552B-CT
Freezer Rack/Box (Quantity)	16/400
Electrical Data	
Power Supply(V/Hz)	110-230V/50Hz-60Hz
Rated Current (A)	9.09A